

PRODUCT DATA SHEET

Bioworld Technology,Inc.

GNB2 polyclonal antibody

Catalog: BS71216 Host: Rabbit Reactivity: Human, Mouse, Rat

BackGround:

Heterotrimeric guanine nucleotide-binding proteins (G proteins), which integrate signals between receptors and effector proteins, are composed of an alpha, a beta, and a gamma subunit. These subunits are encoded by families of related genes. This gene encodes a beta subunit. Beta subunits are important regulators of alpha subunits, as well as of certain signal transduction receptors and effectors. This gene contains a trinucleotide (CCG) repeat length polymorphism in its 5' UTR.

Product:

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

37kDa

Swiss-Prot:

P62879

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB,1:500 - 1:2000|IHC,1:50 - 1:200|IF/ICC,1:50 - 1:200


Storage&Stability:

Store at $4 \, \mathbb{C}$ short term. Aliquot and store at $-20 \, \mathbb{C}$ long term. Avoid freeze-thaw cycles.

Category:

Polyclonal Antibodies

DATA:

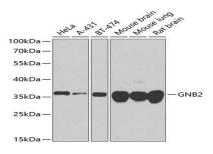

Western blot analysis of extracts of various cell lines, using GNB2 antibody at 1:1000 dilution.

Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.

br/>Lysates/proteins: 25ug per lane.

lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit .

Exposure time: 90s.

Immunohistochemistry of paraffin-embedded rat stomach using GNB2 Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Immunohistochemistry of paraffin-embedded mouse testis using GNB2 Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Immunohistochemistry of paraffin-embedded mouse pancreas using GNB2 Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u> Tel: 0086-025-68037686 Fax: 0086-025-68035151