

ERCC5 polyclonal antibody

Catalog: BS6595

Host: R

Rabbit

Reactivity: Human, Mouse, Rat

BackGround:

This gene encodes a single-strand specific DNA endonuclease that makes the 3' incision in DNA excision repair following UV-induced damage. The protein may also function in other cellular processes, including RNA polymerase II transcription, and transcription-coupled DNA repair. Mutations in this gene cause xeroderma pigmentosum complementation group G (XP-G), which is also referred to as xeroderma pigmentosum VII (XP7), a skin disorder characterized by hypersensitivity to UV light and increased susceptibility for skin cancer development following UV exposure. Some patients also develop Cockayne syndrome, which is characterized by severe growth defects, mental retardation, and cachexia. Read-through transcription exists between this gene and the neighboring upstream BIVM (basic, immunoglobulin-like variable motif containing) gene.

Product:

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

Refer to figures

Swiss-Prot:

P28715

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

IF/ICC,1:50 - 1:200

Storage&Stability:

Store at $4 \,^{\circ}{\rm C}$ short term. Aliquot and store at $-20 \,^{\circ}{\rm C}$ long term. Avoid freeze-thaw cycles.

Category:

Polyclonal Antibodies

Bioworld Technology, Inc.


 Add:
 1660 South Highway 100, Suite 500 St. Louis Park, MN 55416,USA.

 Email:
 info@bioworlde.com


 Tel:
 6123263284

 Fax:
 6122933841


DATA:

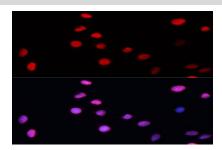

Immunofluorescence analysis of C6 cells using ERCC5 Polyclonal Antibody at dilution of 1:100. Blue: DAPI for nuclear staining.

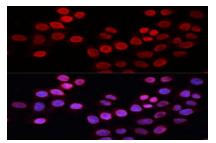

Immunofluorescence analysis of HeLa cells using ERCC5 Polyclonal Antibody at dilution of 1:100. Blue: DAPI for nuclear staining.

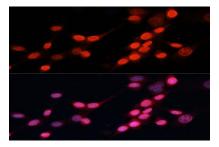

Immunofluorescence analysis of NIH-3T3 cells using ERCC5 Polyclonal Antibody at dilution of 1:100. Blue: DAPI for nuclear staining.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld technology, co. Ltd. Add: No 9, weidi road Qixia District Nanjing, 210046, P. R. China. Email: info@biogot.com Tel: 0086-025-68037686 Fax: 0086-025-68035151