

PRODUCT DATA SHEET

Bioworld Technology,Inc.

Insulin polyclonal antibody

Catalog: BS6890 Host: Rabbit Reactivity: Human, Mouse, Rat

BackGround:

After removal of the precursor signal peptide, proinsulin is post-translationally cleaved into three peptides: the B chain and A chain peptides, which are covalently linked via two disulfide bonds to form insulin, and C-peptide. Binding of insulin to the insulin receptor (INSR) stimulates glucose uptake. A multitude of mutant alleles with phenotypic effects have been identified. There is a read-through gene, INS-IGF2, which overlaps with this gene at the 5' region and with the IGF2 gene at the 3' region. Alternative splicing results in multiple transcript variants.

Product:

1 mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

6KDa

Swiss-Prot:

P01308

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB,1:500 - 1:2000|IHC,1:50 - 1:200|IF/ICC,1:50 - 1:200

Storage&Stability:

Store at $4 \,\mathrm{C}$ short term. Aliquot and store at $-20 \,\mathrm{C}$ long term. Avoid freeze-thaw cycles.

Category:


Polyclonal Antibodies

DATA:


Western blot analysis of various lysates, using Insulin antibody at

1:671 dilution.
Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.
dr/>Lysates/proteins: 25ug per lane.
br/>Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit .
Exposure time: 30s.

Immunohistochemistry of paraffin-embedded rat pancreatic islet using Insulin Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

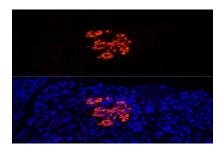

Immunofluorescence analysis of human pancreas cells using Insulin Rabbit pAb at dilution of 1:50 . Blue: DAPI for nuclear staining.

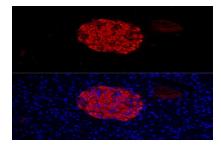

Immunofluorescence analysis of mouse pancreas cells using Insulin Rabbit pAb at dilution of 1:50 . Blue: DAPI for nuclear staining.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841

Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151