

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# VCP polyclonal antibody

Catalog: BS7279 Host: Rabbit Reactivity: Human, Mouse, Rat

#### **BackGround:**

The protein encoded by this gene is a member of a family that includes putative ATP-binding proteins involved in vesicle transport and fusion, 26S proteasome function, and assembly of peroxisomes. This protein, as a structural protein, is associated with clathrin, and heat-shock protein Hsc70, to form a complex. It has been implicated in a number of cellular events that are regulated during mitosis, including homotypic membrane fusion, spindle pole body function, and ubiquitin-dependent protein degradation.

#### **Product:**

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

#### **Molecular Weight:**

90KDa

#### **Swiss-Prot:**

P55072

#### **Purification&Purity:**

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

#### **Applications:**

WB,1:500 - 1:2000|IHC,1:50 - 1:200|IF/ICC,1:50 - 1:200


#### Storage&Stability:

Store at  $4\,\mathrm{C}$  short term. Aliquot and store at -20  $\mathrm{C}$  long term. Avoid freeze-thaw cycles.


#### **Category:**

Polyclonal Antibodies


#### DATA:

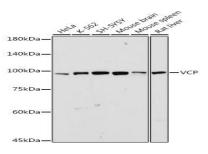

Western blot analysis of extracts of various cell lines, using VCP antibody at 1:1000 dilution.<br/>
Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.<br/>
br/>Lysates/proteins: 25ug per lane.<br/>
lane.<br/>
Blocking buffer: 3% nonfat dry milk in TBST.<br/>
Detection: ECL Basic Kit .<br/>
Exposure time: 30s.

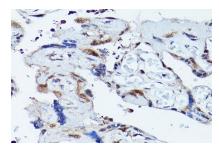

Immunohistochemistry of paraffin-embedded human placenta using VCP Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

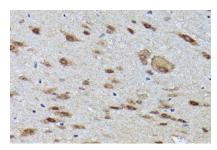

Immunohistochemistry of paraffin-embedded mouse spinal cord using VCP Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

## Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841

## Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,


P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



## PRODUCT DATA SHEET

Bioworld Technology,Inc.

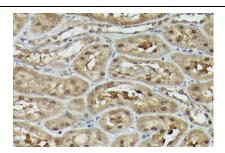

Immunohistochemistry of paraffin-embedded rat kidney using VCP Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

#### Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Tel: <a href="mailto:0086-025-68037686">0086-025-68037686</a>
Fax: <a href="mailto:0086-025-68035151">0086-025-68035151</a>