

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# **VPS18** polyclonal antibody

Catalog: BS76891 Host: Rabbit Reactivity: Human, Mouse, Rat

#### **BackGround:**

Vesicle mediated protein sorting plays an important role in segregation of intracellular molecules into distinct organelles. Genetic studies in yeast have identified more than 40 vacuolar protein sorting (VPS) genes involved in vesicle transport to vacuoles. This gene encodes the human homolog of yeast class C Vps18 protein. The mammalian class C Vps proteins are predominantly associated with late endosomes/lysosomes, and like their yeast counterparts, may mediate vesicle trafficking steps in the endosome/lysosome pathway.

#### **Product:**

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

#### **Molecular Weight:**

110kDa

#### **Swiss-Prot:**

Q9P253

## **Purification&Purity:**

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

## **Applications:**

WB,1:500 - 1:2000|IF/ICC,1:50 - 1:200


#### Storage&Stability:

Store at  $4 \, \mathbb{C}$  short term. Aliquot and store at  $-20 \, \mathbb{C}$  long term. Avoid freeze-thaw cycles.


#### **Modification:**

Unmodification


#### DATA:

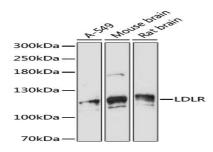

Western blot analysis of extracts of various cell lines, using VPS18 antibody at 1:1000 dilution.<br/>
Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.<br/>
br/>Lysates/proteins: 25ug per lane.<br/>
lane.<br/>
Blocking buffer: 3% nonfat dry milk in TBST.<br/>
Detection: ECL Basic Kit .<br/>
Exposure time: 4min.

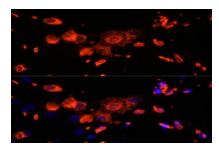

Immunofluorescence analysis of C6 cells using Vps18 antibody at dilution of 1:100. Blue: DAPI for nuclear staining.

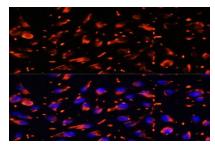

Immunofluorescence analysis of L929 cells using Vps18 antibody at dilution of 1:100. Blue: DAPI for nuclear staining.

## Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,


P. R. China.

Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



## PRODUCT DATA SHEET

Bioworld Technology,Inc.

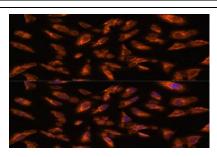

Immunofluorescence analysis of U-2 OS cells using Vps18 antibody at dilution of 1:100. Blue: DAPI for nuclear staining.

### Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151