

## ZAK polyclonal antibody

Catalog: BS8865

Host:

Rabbit

Reactivity: Human, Mouse, Rat

## **BackGround:**

This gene is a member of the MAPKKK family of signal transduction molecules and encodes a protein with an N-terminal kinase catalytic domain, followed by a leucine zipper motif and a sterile-alpha motif (SAM). This magnesium-binding protein forms homodimers and is located in the cytoplasm. The protein mediates gamma radiation signaling leading to cell cycle arrest and activity of this protein plays a role in cell cycle checkpoint regulation in cells. The protein also has pro-apoptotic activity. Alternate transcriptional splice variants, encoding different isoforms, have been characterized.

#### **Product:**

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

**Molecular Weight:** 

55kDa/91kDa/100kDa

**Swiss-Prot:** 

## **Purification&Purity:**

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

#### **Applications:**

WB,1:500 - 1:2000|IHC,1:50 - 1:200|IF/ICC,1:20 - 1:50


## Storage&Stability:

Store at  $4 \,^{\circ}{\rm C}$  short term. Aliquot and store at  $-20 \,^{\circ}{\rm C}$  long term. Avoid freeze-thaw cycles.


**Category:** 

## Polyclonal Antibodies


**DATA:** 

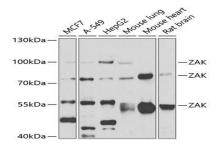

Western blot analysis of extracts of various cell lines, using ZAK antibody at 1:1000 dilution.<br/>Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.<br/>br/>Lysates/proteins: 25ug per lane.<br/>br/>Blocking buffer: 3% nonfat dry milk in TBST.<br/>Detection: ECL Basic Kit .<br/>br/>Exposure time: 90s.

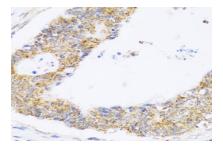

Immunohistochemistry of paraffin-embedded human colon carcinoma using ZAK antibody at dilution of 1:100 .Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

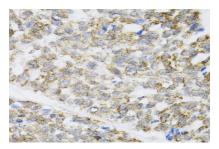

Immunohistochemistry of paraffin-embedded human esophageal cancer using ZAK antibody at dilution of 1:100 .Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

#### Bioworld Technology, Inc.

 Add:
 1660 South Highway 100, Suite 500 St. Louis Park, MN 55416,USA.

 Email:
 info@bioworlde.com

 Tel:
 6123263284


 Fax:
 6122933841

## Bioworld technology, co. Ltd. Add: No 9, weidi road Qixia District Nanjing, 210046, P. R. China. Email: info@biogot.com Tel: 0086-025-68037686 Fax: 0086-025-68035151



# **PRODUCT DATA SHEET**

Bioworld Technology,Inc.

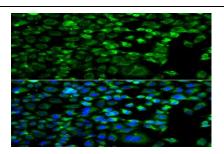

Immunofluorescence analysis of A549 cells using ZAK antibody . Blue: DAPI for nuclear staining.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.		Bioworld technology, co. Ltd.	
Add:	1660 South Highway 100, Suite 500 St. Louis Park,	Add: No 9, weidi road Qixia District Nanjing, 210	046,
	MN 55416,USA.	P. R. China.	
Email:	info@bioworlde.com	Email: <u>info@biogot.com</u>	
Tel:	6123263284	Tel: 0086-025-68037686	
Fax:	6122933841	Fax: 0086-025-68035151	