

Bioworld Technology,Inc.

ILK polyclonal antibody

Catalog: BS90712

Host: R

Rabbit

Reactivity: Human, Mouse, Rat

BackGround:

Integrins are heterodimers composed of non-covalently associated transmembrane a and b subunits. The 16 a and 8 b subunits heterodimerize to produce more than 20 different receptors. Most integrin receptors bind to ligands that are components of the extracellular matrix. Certain integrins can also bind to soluble ligands such as Fibrinogen, or to counterreceptors on adjacent cells, such as the intracellular adhesion molecules (ICAMs), leading to aggregation of cells. In addition to mediating cell adhesion and cytoskeletal organization, integrins function as signaling receptors. Signals transduced by integrins play a role in many biological processes, including cell growth, differentiation, migration and apoptosis. ILK (integrin-linked kinase) was identified as a serine/threonine kinase that phosphorylates b1 and b3 integrins. ILK expression has been shown to be reduced in response to Fibronectin, a known integrin ligand. Overexpression of ILK was shown to upregulate the Fibronectin matrix assembly in epithelial cells, indicating a potential role for ILK in cell growth, cell survival and tumorigenesis.

Product:

Rabbit IgG, 1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

51 kDa

Swiss-Prot:

Q13418(Human) O55222(Mouse) Q99J82(Rat)

Purification&Purity:

ProA affinity purified

Applications:

WB:1:1,000-1:2,000 ICC:1:50-1:200 IHC:1:50-1:200 FC:1:50-1:100


Storage&Stability:

Store at +4 $^{\circ}$ C after thawing. Aliquot store at -20 $^{\circ}$ C or -80 $^{\circ}$ C. Avoid repeated freeze / thaw cycles.

Specificity:


ILK polyclonal antibody detects endogenous levels of ILK protein.

DATA:

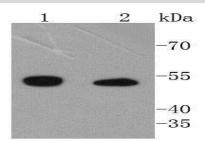

Western blot analysis of Integrin linked ILK on different lysates using anti-Integrin linked ILK antibody at 1/1,000 dilution. Positive control: Lane 1: Hela

Lane 2: 293T

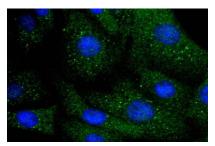

ICC staining Integrin linked ILK in NIH/3T3 cells (green). The nuclear counter stain is DAPI (blue). Cells were fixed in paraformaldehyde, permeabilised with 0.25% Triton X100/PBS.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

 Add:
 1660 South Highway 100, Suite 500 St. Louis Park, MN 55416,USA.

 Email:
 info@bioworlde.com

 Tel:
 6123263284

 Fax:
 6122933841

Bioworld technology, co. Ltd.

 Add:
 No 9, weidi road Qixia District Nanjing, 210046, P. R. China.

 Email:
 info@biogot.com

 Tel:
 0086-025-68037686

 Fax:
 0086-025-68035151