

# PRODUCT DATA SHEET

Bioworld Technology,Inc.

# **Neutrophil Elastase polyclonal antibody**

Catalog: BS90933 Host: Rabbit Reactivity: Human

### **BackGround:**

Neutrophil Elastase (NE) is a serine protease that is expressed in bone marrow precursor cells, stored in peripheral blood granulocytes, and implicated in the progression of a variety of inflammatory diseases, including idiopathic pulmonary fibrosis, rheumatoid arthritis, adult respiratory distress syndrome and cystic fibrosis. In neutrophils, Neutrophil Elastase contributes largely to the proteolysis of phagocytosed proteins, the migration of neutrophils and the remodeling of tissues following injury. Neutrophil Elastase, which is also designated medullasin, is secreted into the extracellular matrix, where it is then capable of destroying connective tissue proteins, including elastin, proteoglycans and Type IV Collagens. Neutrophil Elastase also mediates proteolysis by cleaving proteins that are associated with the complement system, such as antithrombin and Fibrinogen. Additionally, Neutrophil Elastase functions as a potent platelet agonist, where it potentiates the aggregation, secretion and mobilization of calcium in response to cathepsin G binding to platelet surface receptors.

## **Product:**

Rabbit IgG, 1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

### **Molecular Weight:**

29 kDa

#### **Swiss-Prot:**

P08246(Human)

## **Purification&Purity:**

ProA affinity purified

## **Applications:**

WB:1:1,000 ICC:1:50

IHC:1:50-1:200 FC:1:10-1:50


## **Storage&Stability:**

Store at +4  $^{\circ}$ C after thawing. Aliquot store at -20  $^{\circ}$ C or -80  $^{\circ}$ C. Avoid repeated freeze / thaw cycles.


# **Specificity:**

Neutrophil Elastase polyclonal antibody detects endogenous levels of Neutrophil Elastase protein.

#### **DATA:**

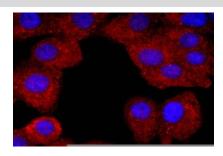

ICC staining Neutrophil Elastase in A549 cells (red). The nuclear counter stain is DAPI (blue). Cells were fixed in paraformaldehyde, permeabilised with 0.25% Triton X100/PBS.

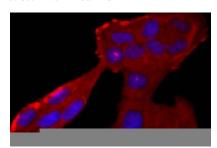

ICC staining Neutrophil Elastase in Hela cells (red). The nuclear counter stain is DAPI (blue). Cells were fixed in paraformaldehyde, permeabilised with 0.25% Triton X100/PBS.

#### Note

For research use only, not for use in diagnostic procedure.

## Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u> Tel: 0086-025-68037686 Fax: 0086-025-68035151